# ReFa

ReFa BEAUTECH SPOT  
Yellow/Blue/Red

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ReFa

ReFa BEAUTTECH SPOT  
Yellow/Blue/Red



## Quick Guide

Household Beauty Device

Model Number: RE-AK16A/RE-AK12A/RE-AK01A

- Make sure to read "Safety Precautions" and "Precautions for Use" in the separate "Instruction Manual" before use.

## Names of the Parts

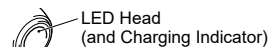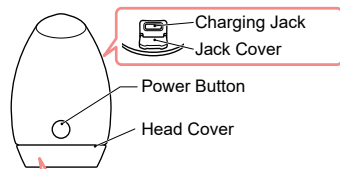When the power is turned on, the two types of LED for skin care light up.  
(1) Near infrared LED (invisible light)  
(2) Visible LED

The emission color varies depending on the product (model number).

- Yellow (RE-AK16A) : Yellow
- Blue (RE-AK12A) : Blue
- Red (RE-AK01A) : Red

- During charging, the low light (low-intensity emission) indicates the charging condition.


\* For details, refer to "Names of the Parts and Functions" in the separate "Instruction Manual".

## Charging

- The AC adapter is not included.
- To charge the main unit, use the included USB charging cable (USB Type-C). An accident and/or fire may result if anything else is used.
- When charging the product, check the rated output of the USB connector (Type-A terminal) connecting destination (OUTPUT) and the rated input of this product (5 V DC, 0.5 A).
- \* Note that a rated output lower than the rated input may cause the connection destination to be raised to a high temperature or may damage it.
- This product operates on a built-in rechargeable battery. When using the product for the first time, charge the battery before using it.
- You may not be able to use the product after a long period of disuse because the remaining charge of the battery is low. If this happens, charge the device again before using it.

- Charge in a location with a temperature range of 10°C to 35°C. Battery performance may deteriorate at lower or higher temperatures due to the nature of the rechargeable battery.
- Cannot be used while charging.


## 1 Open the jack cover downward and insert the included USB charging cable's Type-C terminal into the charging jack of the main unit

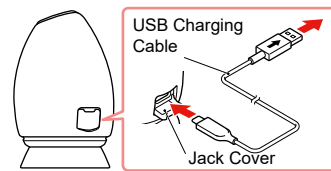- When charging starts, a short beep will sound.
- When charging, place the head cover upside down to use as a stand, and place the product on the upside-down head cover to make it easy to check the lighting of the LED head (which also functions as the charging indicator).

## 2 Check the LED head for the charging condition

- Charging: Slow blinking (low light)
- Charging complete (fully charged): Not lit

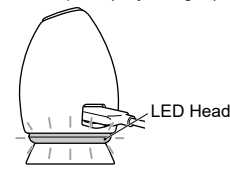- When charging is complete, short beeps will also sound.
- The charging time varies by the remaining charge of the battery.
- The available number of uses that the fully charged battery allows is as follows.
  - Yellow (RE-AK16A): Approx. 16
  - Blue (RE-AK12A): Approx. 10
  - Red (RE-AK01A): Approx. 10
- When the remaining charge of the battery becomes low during use, short beeps will sound.
- When the remaining charge of the battery runs out during use, a long beep will sound and then the power will go off.

## Usage Area

The recommended usage area is the face.

\* You can also use this product on your body.

## Basic Usage


- This product is not waterproof.

## 1 Remove your makeup and wash your face

- We recommend you use the product after you have washed your face or after applying cosmetic water.

## 2 Remove the head cover

## 3 Press and hold the power button

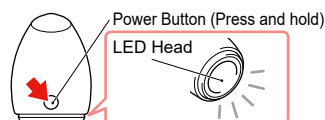- The LED head will light up and a melody will also sound.

## 4 Press the power button to select the vibration mode according to your preference

- A short beep sounds and the status changes between "vibration off" and "vibration on" alternately each time the power button is pressed.
- Note that the power goes off if you press and hold the power button.

## 5 Apply the bottom surface of the LED head so that it lightly touches the skin

- Apply the LED light to areas of concern on your skin for a maximum of about 5 minutes per area.
- When about 5 minutes have elapsed, a long beep will sound to report the approximate time.

- Do not use it on the eyelids, beneath the brows, eyeballs, and the area near the eyeballs.
- When applying the LED head to the areas where the LED head can point toward the eyes, such as the cheeks and the forehead, close your eyes to prevent from directly looking into the LED light.
- Do not press hard when applying the LED head to the skin. It may put stress on your skin.

## 6 Press and hold the power button to end care

- This product has an automatic power-off function. The power will automatically go off when approximately 10 minutes have elapsed after power-on, even if you do not press and hold the power button.

## 7 Wipe off any dirt or grime on the surface of the LED head and attach the head cover

- Wipe it softly with a tissue or a cloth dampened with water that has been thoroughly wrung out.

- Do not wash the main unit with water.
- Turn the power off when performing maintenance after use.
- Do not use paint thinner, benzene, alcohol, nail polish remover, or other such materials on the product. There is a risk of deformation or cracking of parts and lead to product damage.
- When the product is not in use, make sure to attach the head cover.

## When the Following Happens

## Q Can I use the product with my makeup on?

- A. Some cosmetics block or reduce the light, so we recommend you use the product after you have washed your face or after applying cosmetic water.

## Q Can I use this product every day?

- A. You can use it every day, but you should not apply it to the same area more than twice a day (about 5 minutes each time).

## Q The LED head warms up while the product is in use.

- A. This is because electronic parts generate heat. It is not a failure, so use the product with peace of mind.

## For Other Product Information

For details on handling the product, such as "Maintenance," "Q&A (frequently asked questions)," or "Safety Precautions," etc., refer to the separate "Instruction Manual."

Publication Date: 2021/11 1st Edition

ReFa

ReFa BEAUTTECH SPOT  
Yellow/Blue/Red

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