ReFa

ReFa CAXA

GUIDE BOOK

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ENGLISH

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
Thank you for purchasing ReFa CAXA.

ReFa represents a sincere desire to see women harness the beauty that flows from within and without them——.

It is a brand that expresses a wish to see you delight in discovering your true radiance with just a minimum of effort.

It also symbolizes an intention to produce a spontaneous smile on your face as you enjoy a sense of fulfillment from simply using our product.

ReFa is here for your everyday pleasure.

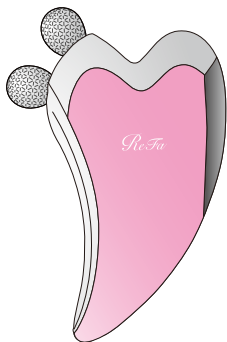■ Esthetic Method

Fascia are targeted through a scooping-up action that reproduces the "friction" technique of estheticians along with a grip-and-glide action that reproduces their "kneading" technique.

■ Microcurrent

The solar panel takes in light and releases microcurrent.


■ Platinum Coating

To protect your skin during contact, the surface of the device's rollers are coated with platinum. This has the added benefit of providing the device with a lasting shine that is pleasing to the eye.

■ Waterproof Construction

Conforms with JIS standards (equivalent to IPX7). So the rollers can be used in the bathtub.

Names of the Parts

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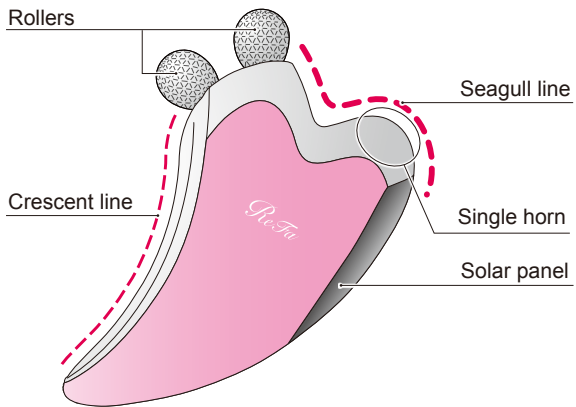ENGLISH

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
■ Important Points in Using the Crescent and Seagull Lines

While making firm contact with your skin, move as if scooping upward.

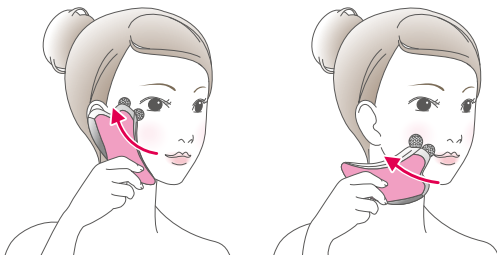
■ Important Points in Using the Rollers

The strength with which the roller pulls up your skin can be adjusted by changing the angle at which you hold the body.

Gradually lower the body and assess the lifting force of the rollers.

*Raise the body if you begin to feel pain.

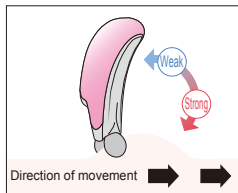How To Use


■ Refer to the instructions in the following pages and use according to your preference.

Please be sure to read the instruction manual before use


- When using this device for the first time, roll it gently a few times over an area of the body that is not very noticeable such as your upper arm, and make sure that there are no abnormal skin reactions after some time has elapsed.
- Persons with metal allergies, or who have conditions such as skin infections or dermatitis, should not use the product.
- Women in the early unstable stages of pregnancy or immediately after childbirth, or women who are pregnant and not in good health should not use this product.
- Remove any accessories (necklaces, piercings, etc.) during use, as they may get caught in the rollers.
- Tie your hair up or wear a hair band when using the device so your hair does not get caught.
- Do not use along the spine or in other locations where bones protrude.
- Please do not use anywhere other than as specified in this manual.

CAXA Work


- (1) Push to give care along the inside of your collarbone with the single horn. Press inside to outside at 3 points each on the left & right sides.

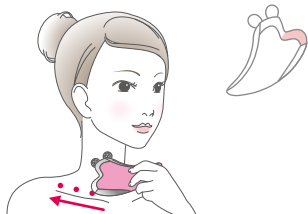
- (2) Fit the seagull line onto your neck muscle, and push and drain toward your collarbone.

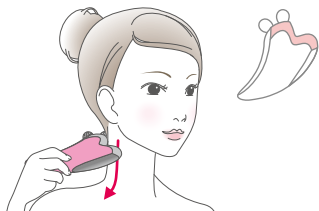
- (3) Fit the seagull line right on your jaw line and scoop up along the face line. Exert force as you do this.

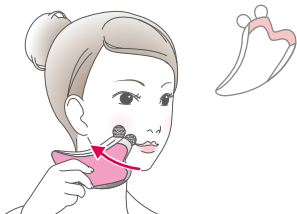
- (4) Place the seagull line below the cheekbone and scoop up from the side of the nose to the front of the ear.

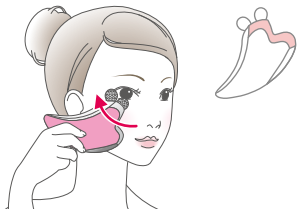- (5) Apply the crescent line tightly to your entire cheek, and scoop up towards your ear. Imagine it lifting your whole cheek.

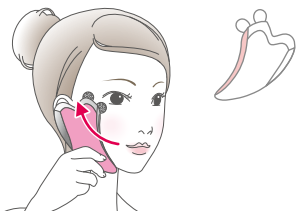
- (6) Apply the crescent line on half of your forehead and pull up diagonally.

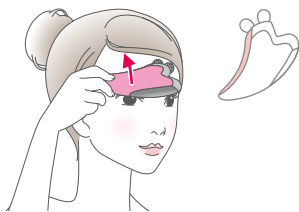
*Make sure your hair does not get caught in the product.

Rolling Care


- (1) Roll left and right with your collarbone between the rollers.

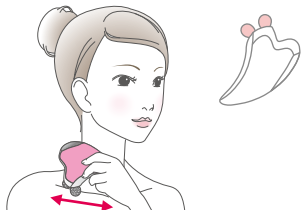
- (2) Roll up and down along the side of the neck.

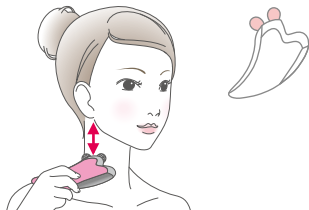
(3) Roll along the face line from the tip of the chin to the bottom of the ear.

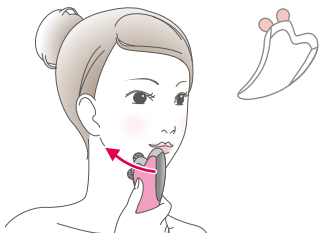(4) Roll the rollers gently along the side of your nose, from the wing of your nose to the inner corner of your eye.

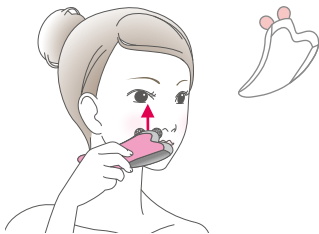(5) Roll the rollers up and down around your eyes. Roll in several spots around your eyes.

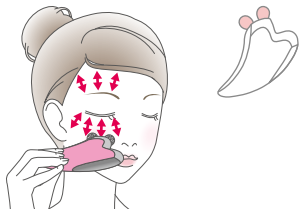(6) Roll the rollers up and down your forehead. Use in a zigzag motion over your entire forehead.

*Make sure your hair does not get caught in the product.

ReFa

ReFa CAXA
GUIDE BOOK
Simplified Chinese

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